


HANAROSTENT®
Biliary Flap Lasso (CCC)

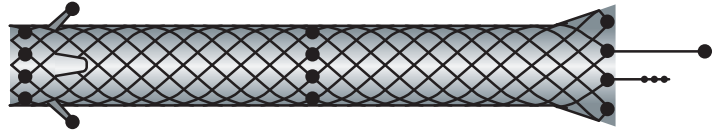

Indication :

- Palliative treatment of bile duct malignant and/or benign strictures

Features :

- Flare on the proximal and 4 Anchoring Flaps for Anti-Migration
- Fully covered stent for Endoscopic Application
- Radiopaque markers for fluoroscopic guidance
- Lasso on the proximal for stent removal and long lasso for Inside-Out stent removal
- Endoscopic Application

Endoscopic Application

	Stent (mm)			Delivery device	
	Diameter	Usable length	Total length*	Length(mm)	Diameter(mm/fr)
BCT-08-040-180	8-11.5	23	40	1800	3.4/10.2
BCT-08-120-180		103	120	1800	3.4/10.2
BCT-10-040-180	10-13.5	23	40	1800	3.4/10.2
BCT-10-120-180		103	120	1800	3.4/10.2

*Increment of 10mm

	(NNN) Non covered	<u>XXX-00-000-000</u>	<u>XXX 00 000-X000</u>
	(CCC) Fully covered	① ② ③ ④	① ② ③ ④
	(NCN) Partially covered	① Product Name	③ Stent Length
		② Stent Diameter	④ Delivery Device Length